

RJN系列软胶囊机生产线

RJN Series Softgel Encapsulation Machine Line

► 用途 (Application)

软胶囊机可将油类液体、混悬液或糊状物定量注入胶膜中，形成大小、形状、色彩各异的密封软胶囊。由于软胶囊具有崩解速度快、生物利用率高、外形雅致、易于吞服、不易氧化、便于贮存和携带等优点，越来越受到消费者的欢迎。软胶囊机由于生产率高、装量准确、质量稳定、成品率高而越来越受到厂家的青睐，现已广泛用于医药、保健品、化妆品和游戏彩弹等领域。


软胶封装系统由主封装机以及输送装置、熔胶罐、配料罐、冷水机、定型转笼干燥机、电气控制面板、带小车的烘箱托盘、储罐等组成。

Softgel encapsulation machine can quantitatively inject the oil, mixing suspension and pasty mass into the capsule to form all kinds of softgel encapsulation with different size, shape and color. Because of the merits of high speed disintegration, easy to oral, hard to oxygenation, easy to store and take, the softgel encapsulation is welcomed by more and more customers. Owing to the high productivity, charging exactly, quality stabilization, high finished product rate, it is favoured by more and more manufacturers and used in the fields of medicine, health products, cosmetic and game-color ball and so on.

The softgel encapsulating system consists of main encapsulation machine as a crucial apparatus and conveyer, gelatin melting tank, dosage tank, water chilling unit, tumble dryer for formalization, electrical control panel, oven tray with dolly, storage tank, etc.



软胶丸工艺流程 Workflow of Soft Gelatin Capsule

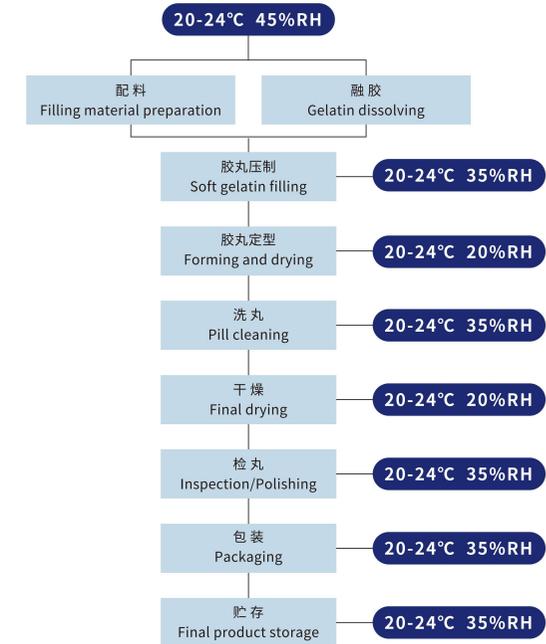

► 技术参数 (Technical Parameters)

参数ITEM	型号TYPE	RJN-120	RJN-200	RJN-250	RJN-300	RJN-320
滚模转速 Rotary Speed of Die Roll	rpm	0-5	0-5	0-5	0-5	0-5
		变频无级调速 Frequency conversion stepless speed regulation				
最大产量 Max. Output	pcs/h	18000	65000	92000	110000	150000
流体单柱塞供给量 Supplying Volume for Flow Single-pillar Pisto	g	0-2	0-2	0-2	0-2	0-2
装量差异 Difference Among Filling Volume	%	±2	±2	±2	±2	±2
整机功率 Power for Whole Machine	kw	4	7.5	9.5	17.6	20.6

